


Dental Handpiece Maintenance Equipment



# Lubrina 2

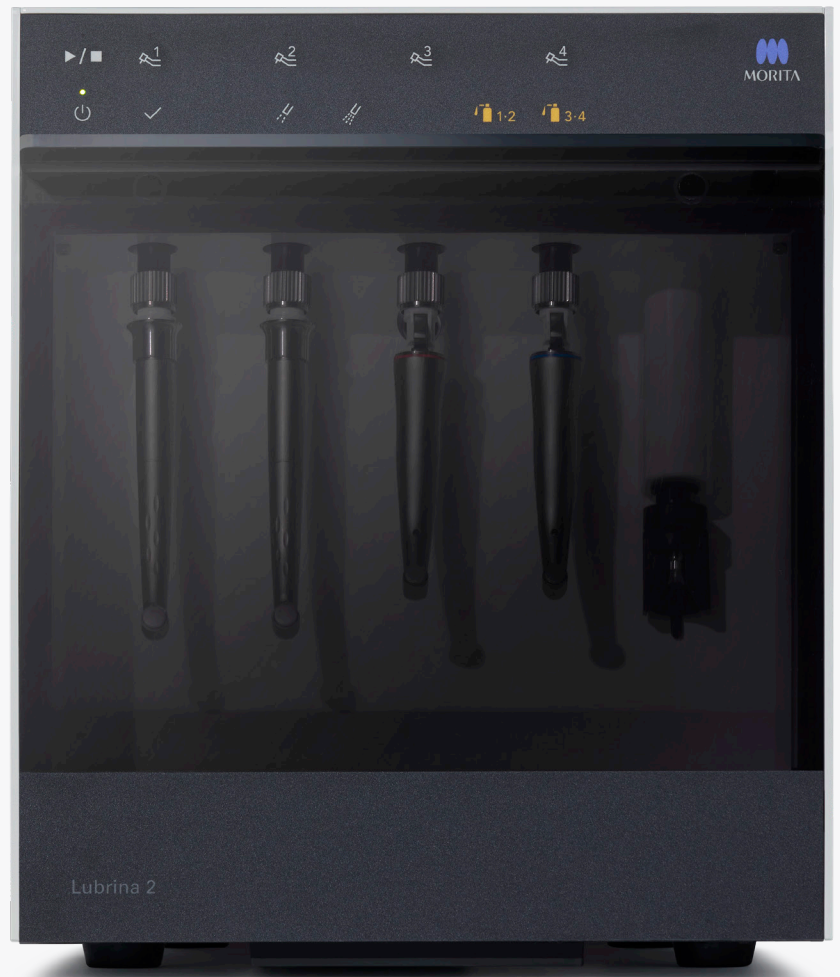

Thinking ahead. Focused on life.

# Seamless Auto-Maintenance

Fast and clean operation, plus compatibility with multiple brands of spray lubricant

Lubrina 2 automatically performs maintenance and lubrication of handpieces which are essential for daily practice. Because the lubrication amount and service time suitable for each handpiece are programmed into Lubrina 2, it can automatically perform maintenance of up to four handpieces in a short time. The clean air-blow system implements lubrication and purging of the inside of the handpiece thoroughly, providing consistent, quality maintenance that doesn't rely on the experience of staff which can occur during manual lubrication. Lubrina 2 accepts multiple brands of spray lubricant.




Approximately 20 seconds per turbine

Oil usage amount is reduced to 1/4



Faster maintenance with excellent cost savings

Since the air-blow system can remove residual water completely using clean air, the efficiency of lubrication is enhanced with oil usage only 1/4 of that compared to manual lubrication. Daily maintenance cost is therefore reduced. In addition, the system can remove surplus oil effectively, providing effective maintenance in mere seconds.\*1

Greatly shortens the maintenance process

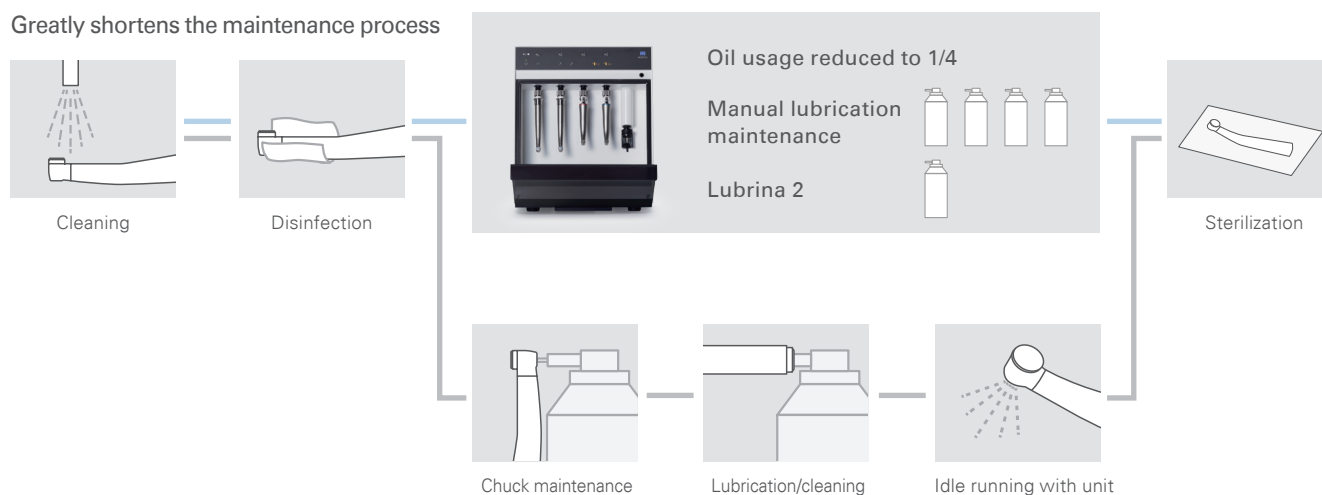

Maintenance of one turbine in approximately 20 seconds

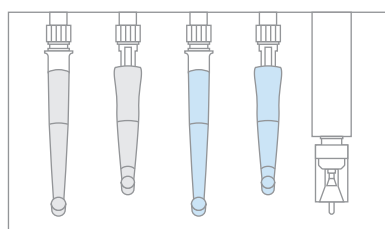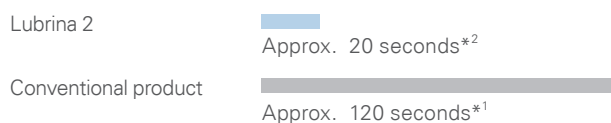

Example: One turbine and contra angle by Morita and one turbine and contra angle by another manufacturer.



### Handpiece and spray compatibility

By replacing the coupling, Lubrina 2 facilitates the maintenance of Tri Auto ZX2 and Dentaport ZX contra-angles, air driven handpieces (ISO 9168/Type 3 joint) as well as all contra-angles and straight handpieces (ISO 3964/E-Type joint). Lubrina 2 is compatible with 4-hole, 5-hole, and 6-pin connections. With its dual lubrication system, one or two different brands of spray may be used in the same cycle. Several brands of spray are compatible by using an optional spray holder adapter.

\*1 Comparison with the long mode of our previous Lubrina


\*2 Approximately 40 seconds for contra-angle handpiece

# Clean Air-Blow System


Double-conduit design achieves clean air

The double-conduit design allows the clean air-blow system to achieve clean air. By employing the dual solenoid valves and clean air dedicated conduit, the system can supply clean air without a mix of oil. During the water removal process before lubrication, the system removes the water in the head which maximizes lubrication results. During idling after lubrication, the surplus oil can be effectively removed by clean air.


Lubrina 2 (Clean Air-Blow System)

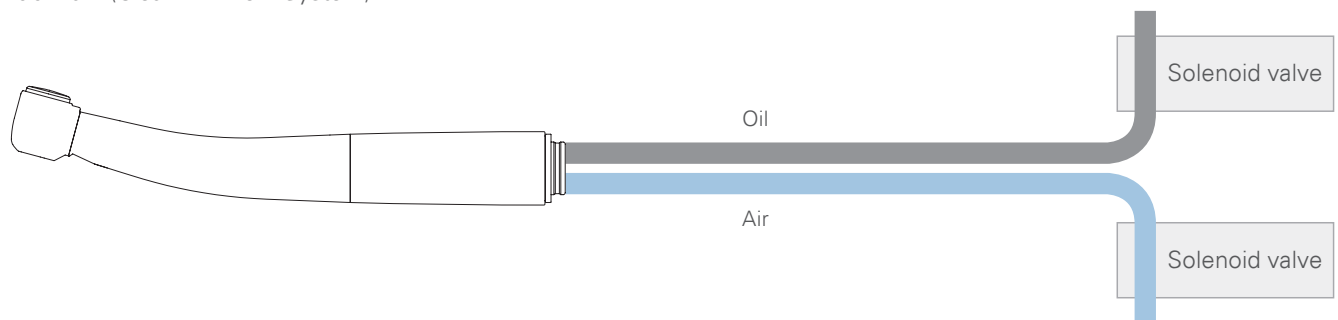

Conventional product\*

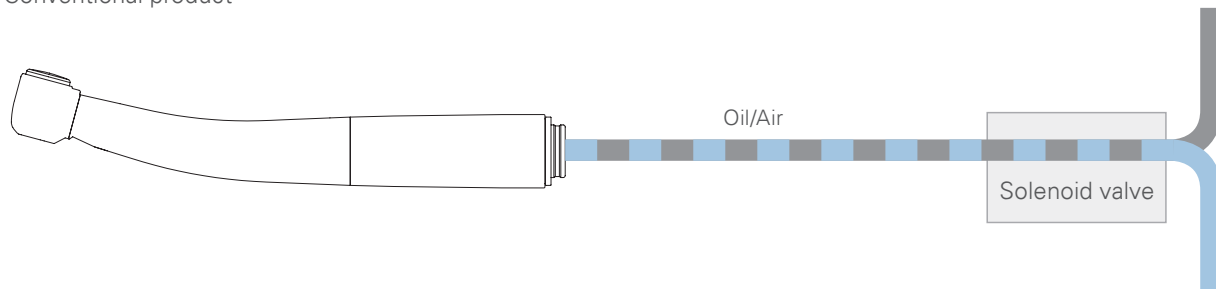


\* Original Lubrina Handpiece Maintenance System

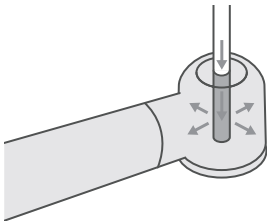


### Easy chuck maintenance

For the turbine and 5x speed contra-angle handpiece, it is essential to lubricate and clean the chuck which holds the bur. The built-in one-touch chuck maintenance mechanism allows for easy maintenance.

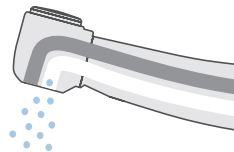


### Water removal during autoclaving

If autoclave is performed while a large amount of water is left in the conduit of the handpiece, steam will not infiltrate and could cause a sterilization defect. To prevent this, clean air (without a mix of oil) is fed into the water feed, chip air conduit to remove the residual water, thereby reducing the possibility of a problem during autoclaving.

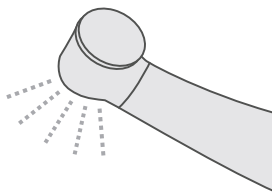


### Surplus oil removal reduces idling time before practice

During handpiece rotation, clean air removes the surplus oil which causes blow-out, and thus reduces the idling time necessary before use.

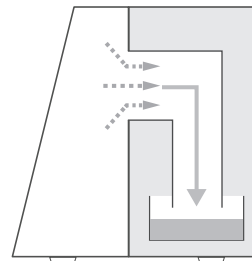


### Reduces scattering of oil mist

The built-in oil mist suction mechanism draws in the oil mist generated during maintenance, thereby preventing the mist from scattering to nearby surfaces and devices. In addition, since oil mist is collected into the oil gathering tray, maintenance of the equipment is also easy.

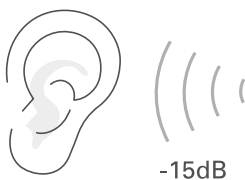

### Reduced operating sound

Lubrina 2 offers a noise reduction of 15dB during operation compared to the original model.

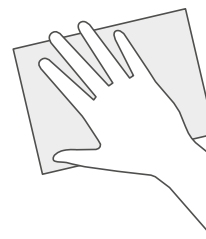


### Easy maintenance

The control panel is flat for easy cleaning, and the front door can also be removed for improved access.

# Features

Compatible with multiple brands

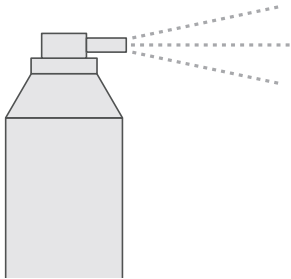

## Oil spray can delivery system









The low viscosity oil mist penetrates into the inside of the handpiece head and delivers effective lubrication.



## Morita Multi Spray (sold separately)

This dedicated oil offers excellent lubrication properties and enhances the performance of TwinPower Turbine and TorqTech.

## Optional Items

Coupling	Target device		
Alpha Coupling Code No.: 7376766 	TwinPower Turbine p, TwinPower Turbine s	4H Coupling Code No.: 7376715 	TwinPower Turbine 4H
KaVo Coupling Code No.: 7376723 	KaVo MULTIflex, KaVo MULTIflex LUX	Sirona Coupling Code No.: 7376731 	Sirona R/F
W&H Coupling Code No.: 7376758 	W&H Roto Quick	LS Coupling Code No.: 7376774 	Handpieces, TorqTech and devices with a connection type that comply with ISO3964 (JIS T5904)
TR-ZX2/DP Coupling Code No.: 7376782 	Tri-Auto ZX2, Tri-Auto mini, Dentaport ZX Root ZX® II	Tri-Auto ZX Coupling Code No.: 7373325 	Tri-Auto ZX

Spray Holder	Matching Spray
Spray Holder (standard)	Morita Multi Spray (LS Spray, AR Spray)
Spray Holder (KaVo/NSK for Pana/Yoshida)	KaVo Spray, NSK PANA SPRAY Plus, YOSHIDA Spray
Spray Holder (for W&H)	W&H SERVICE Oil F1
Spray Holder (for Sirona)	Sirona T1 Spray

Spray	
Morita Multi Spray	

\* The company names and product names are trademarks or registered trademarks of each company.

# Specifications

<b>Brand name:</b>	Lubrina 2
<b>Model:</b>	HIM-2
<b>Power voltage:</b>	AC 100V-240V 50-60Hz
<b>Power consumption:</b>	10-25 VA
<b>Maximum input air pressure:</b>	0.7 MPa max.
<b>Air pressure used:</b>	0.3-0.5 MPa (recommended 0.4 MPa)
<b>Consumed air flow rate:</b>	60-100 NL / min (0.3-0.5 MPa)
<b>Weight:</b>	Approx. 8 kg
<b>External dimensions:</b>	W300×D300×H365mm
<b>Manufactured and sold by:</b>	J. MORITA MFG. CORP.


## Functions

Clean Air-blow System (Dual solenoid valves)  
 Mode 1 (for turbine)  
 Mode 2 (for contra-angle handpiece)  
 Flushing mode  
 Air-blow mode  
 Chuck maintenance mechanism  
 Spray can system  
 Dual spray system  
 Residual water removal mechanism  
 Oil mist suction mechanism  
 Compatible with handpieces manufactured by various companies

Oil absorption pad  
 Oil absorption sheet  
 Oil absorption sheet for front door  
 Sheet stopper for door  
 Power supply cord  
 Air tube  
 Spray Holders (standard) x 2  
 (for Morita Multi Spray)  
 \* Lubrina 2 is not provided with couplings.

## Accessories

Dimensional outline drawing / Installation clearances (mm)

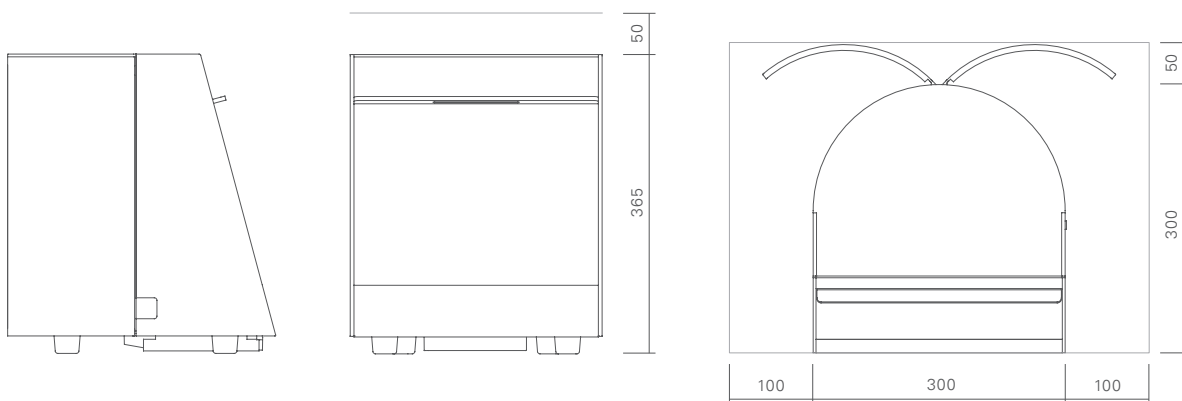



Development and Manufacturing

**J. MORITA MFG. CORP.**

680 Higashihama Minami-cho, Fushimi-ku,  
Kyoto 612-8533, Japan  
T +81. (0)75. 611 2141, F +81. (0)75. 622 4595

**Morita Global Website**

**[www.morita.com](http://www.morita.com)**

Distribution

**J. MORITA CORP.**

3-33-18 Tarumi-cho, Suita-shi, Osaka 564-8650, Japan  
T +81. (0)6. 6380 1521, F +81. (0)6. 6380 0585

**J. MORITA USA, INC.**

9 Mason, Irvine CA 92618, USA  
T +1. 949. 581 9600, F +1. 949. 581 8811

**J. MORITA EUROPE GMBH**

Justus-von-Liebig-Strasse 27b, 63128 Dietzenbach, Germany  
T +49. (0)6074. 836 0, F +49. (0)6074. 836 299

**MORITA DENTAL ASIA PTE. LTD.**

150 Kampong Ampat  
#06-01A KA Centre, Singapore 368324  
T +65. 6779. 4795, F +65. 6777. 2279

**J. MORITA CORP. AUSTRALIA & NEW ZEALAND**

Suite 2.05, 247 Coward Street, Mascot NSW 2020, Australia  
T +61. (0)2. 9667 3555, F +61. (0)2. 9667 3577

**J. MORITA CORP. MIDDLE EAST**

4 Tag Al Roasaa, Apartment 902, Saba Pacha 21311 Alexandria, Egypt  
T +20. (0)3. 58 222 94, F +20. (0)3. 58 222 96

**J. MORITA CORP. INDIA**

Filix Office No.908, L.B.S. Marg, Opp. Asian Paints, Bhandup (West), Mumbai 400078, India  
T +91-22-2595-3482

**J. MORITA MFG. CORP. INDONESIA**

28F, DBS Bank Tower, Jl. Prof. Dr. Satrio Kav. 3-5, Jakarta 12940, Indonesia  
T +62-21-2988-8332, F + 62-21-2988-8201

**SIAMDENT CO., LTD.**

71/10 Mu 5, Thakham, Bangpakong, Chachuengsao 24130, Thailand  
T +66. 38. 573042, F +66. 38. 573043  
[www.siamdent.com](http://www.siamdent.com)

Diagnostic and Imaging Equipment

Treatment Units

Handpieces and Instruments

Endodontic System

Laser Equipment

Laboratory Devices

Educational and Training Systems

Auxiliaries